



ISO 9001:2008 & 13485:2012



# LED SURGICAL LIGHTS

KAMBRIDGE | HEXA | HEXA-CLR-TA



## Kay & Company

[www.kaycoindia.com](http://www.kaycoindia.com)



# Reliable and Effective

The design of the lighthead shape fulfills the requirement of modern laminar flow OR




The open design of the light reduces the surface area that impedes the flow of air from the air supply ceilings, resulting in a less turbulent flow underneath the light.

- The totally enclosed lighthead surface is very easy to decontaminate
- The lenses plate of the lighthead is durable, ageing-resistant and easy to clean
- The handles of lightheads can be easily removed for soaking and autoclavable disinfection


## ✱ Supramaximal effective light field

With high light intensity and large spot size, 1200mm depth of field provides good lighting without having to refocus when the surgical wound location changes. The light head can be positioned for away from the wound to ensure better laminar flow.

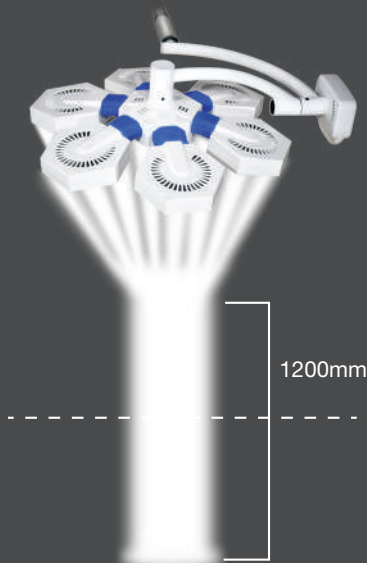


## ✱ Focusing System

The brightness of different LEDs is adjusted automatically to change the spot size, which will avoid the possible occurrence of mechanical malfunctions.

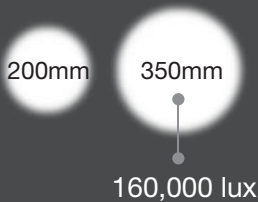


## Endo Mode

To meet the requirement of low ambient illumination.

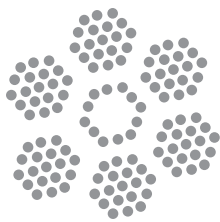

General Mode

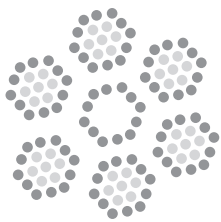

Endo Mode

## Perfect Shadow-free illuminance effect

Thanks to large amount of LED lights, and multiple emitting positions of LED light sources, any obstructions with surgeon's head, shoulder or hands will not affect the surgical wound area with any cast shadows.

### ✳ Control Panel

Touch control panel on the light head.

- On-Off
- Illumination intensity
- Color temperature
- Endo illumination mode




### ✳ Safety Glass for Lens Protection

We offer LED lights with protective fire resistant ESG safety glass cover for scratch proof, fast disinfection process and improve the light penetration. It is made of monolayer safety glass which protects the lens from ultra-violet and infrared rays.



### ✳ Sterilizable Handle

We provide high quality Sterilizable Handle.

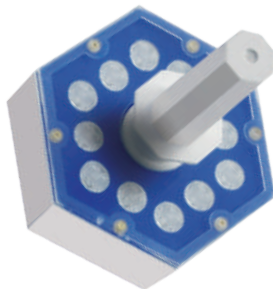

### ✳ HD Documentation for Continuous Medical Education (

- Manual focus
- Integrated central camera or independent third-arm camera (Optional)



## Maximum economy and highest reliability

Its technological advances make it a good long term investment due to low maintenance and operating running cost. With various light head configurations, it maximizes the comfort of the surgeon and forbids the tissue dehydration even during long time operations.

- The minimal life expectancy of each LED is 50,000 hours
- Sophisticated LED thermal management keeps the temperature-rise of the space above surgeon's head and surgical region extremely low.



# INDEX

## KAMBRIDGE LED SERIES

KAMBRIDGE Ceiling LED Lights - Double Light Head	1-5
KAMBRIDGE Ceiling LED Lights - Single Light Head	6-7
KAMBRIDGE Ceiling LED Lights - Mobile Model	6-7

## HEXA LED SERIES

HEXA Ceiling LED Lights - Double Light Head	8-12
HEXA Ceiling LED Lights - Single Light Head	13-14
HEXA Ceiling LED Lights - Mobile Model	13-14

## HEXA-CLR-TA LED SERIES

HEXA-CLR-TA Ceiling LED Lights - Double Light Head	15-17
HEXA-CLR-TA LED Lights - Single Light Head	18-19
HEXA-CLR-TA LED Lights - Mobile Model	18-19

## CLR-TA LED SERIES

CLR-TA Ceiling LED Lights - Double Light Head	20-22
CLR-TA Ceiling LED Lights - Single Light Head	23-24
CLR-TA Ceiling LED Lights - Mobile Model	23-24

## OPERATION & EXAMINATION LIGHTS

AL-HEXA Examination Light	25
AL-C Examination Light	25

NOTE: Reproduction of pictures or specifications in part or full will be considered as illegal and Kay & Company reserves the right to initiate legal action.



KAMBRIDGE 6+6

NUMBER OF LED  
126 + 126

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,60,000 + 1,60,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.



**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 6+5

NUMBER OF LED  
126 + 105

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,60,000 + 1,40,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.



**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 6+4

NUMBER OF LED  
126 + 84

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,60,000 + 1,30,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 6+3

NUMBER OF LED  
126 + 63

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,60,000 + 1,00,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 5+5



NUMBER OF LED  
105 + 105

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,40,000 + 1,40,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 5+4



NUMBER OF LED  
105 + 84

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,40,000 + 1,30,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional



KAMBRIDGE 5+3

**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



**NUMBER OF LED**  
105 + 63

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
1,40,000 + 1,00,000 Lux ±10%



Memory Function



ENDO Mode



Adjustable Colour Temp.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 4+4

**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



**NUMBER OF LED**  
84 + 84

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
1,30,000 + 1,30,000 Lux ±10%



Memory Function



ENDO Mode



Adjustable Colour Temp.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 4+3

NUMBER OF LED  
84 + 63

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,30,000 + 1,00,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 3+3

NUMBER OF LED  
63 + 63

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,00,000 + 1,00,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 6



NUMBER OF LED  
126

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,60,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

MOBILE LED LIGHT

Battery Backup (Optional)



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 5



NUMBER OF LED  
105

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,40,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

MOBILE LED LIGHT

Battery Backup (Optional)



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 4

**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

**NUMBER OF LED**  
84

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
1,30,000 Lux ±10%



Memory Function



ENDO Mode



Adjustable Colour Temp.

**MOBILE LED LIGHT**  
Battery Backup (Optional)



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

KAMBRIDGE 3

**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

**NUMBER OF LED**  
63

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
1,00,000 Lux ±10%



Memory Function



ENDO Mode



Adjustable Colour Temp.

**MOBILE LED LIGHT**  
Battery Backup (Optional)



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	600mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 126+126

NUMBER OF LED  
126 + 126

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,60,000 + 1,60,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 126+105

NUMBER OF LED  
126 + 105

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,50,000 + 1,60,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional



HEXA 126+84

NUMBER OF LED  
126 + 84

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,30,000 + 1,60,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 126+63

NUMBER OF LED  
126 + 63

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,60,000 + 90,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 105+105

NUMBER OF LED  
105 + 105

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,50,000 + 1,50,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 105+84

NUMBER OF LED  
105 + 84

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,50,000 + 1,30,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 105+63

NUMBER OF LED  
105 + 63

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,50,000 + 90,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 84+84

NUMBER OF LED  
84 + 84

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,30,000 + 1,30,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 84+63



NUMBER OF LED  
84 + 63

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,30,000 + 90,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 63+63



NUMBER OF LED  
63 + 63

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
90,000 + 90,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 126

**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

**NUMBER OF LED**  
126

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
1,60,000 Lux ±10%



Memory Function



ENDO Mode



Adjustable Colour Temp.

**MOBILE LED LIGHT**



**MOBILE LED LIGHT**  
with battery backup



Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 105

**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

**NUMBER OF LED**  
105

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
1,50,000 Lux ±10%



Memory Function



ENDO Mode



Adjustable Colour Temp.

**MOBILE LED LIGHT**



**MOBILE LED LIGHT**  
with battery backup



Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional



HEXA 84



**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

**MOBILE LED LIGHT**      **MOBILE LED LIGHT**  
with battery backup



**NUMBER OF LED**  
84

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
1,30,000 Lux ±10%



Memory Function



ENDO Mode



Adjustable Colour Temp.

Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA 63



**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

**MOBILE LED LIGHT**      **MOBILE LED LIGHT**  
with battery backup



**NUMBER OF LED**  
63

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
90,000 Lux ±10%



Memory Function



ENDO Mode



Adjustable Colour Temp.

Easily detachable connector based OT Light / Optional

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA CLR-TA 5+5

NUMBER OF LED  
95 + 95

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,40,000 + 1,40,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA CLR-TA 5+4

NUMBER OF LED  
95 + 76

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,40,000 + 1,10,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA CLR-TA 5+3

NUMBER OF LED  
57 + 95

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
80,000 + 1,40,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA CLR-TA 4+4

NUMBER OF LED  
76 + 76

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,10,000 + 1,10,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA CLR-TA 4+3

NUMBER OF LED  
57 + 76

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
80,000 + 1,10,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA CLR-TA 3+3

NUMBER OF LED  
57 + 57

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
80,000 + 80,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.


SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA CLR-TA 5

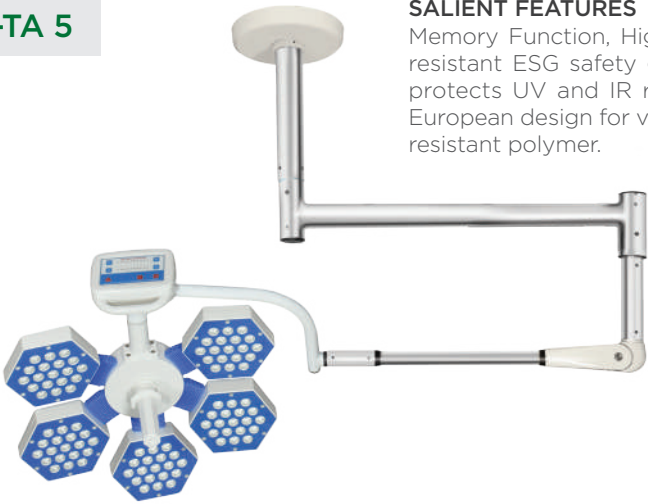

**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

**NUMBER OF LED**  
95

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
1,40,000 Lux ±10%



Memory Function



ENDO Mode



Adjustable Colour Temp.

MOBILE LED LIGHT



MOBILE LED LIGHT with battery backup



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

HEXA CLR-TA 4



**SALIENT FEATURES**  
Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

**NUMBER OF LED**  
76

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
1,10,000 Lux ±10%



Memory Function



ENDO Mode



Adjustable Colour Temp.

MOBILE LED LIGHT



MOBILE LED LIGHT with battery backup



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional



HEXA CLR-TA 3

SALIENT FEATURES

Memory Function, High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

NUMBER OF LED  
57

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
80,000 Lux ±10%



Memory  
Function



ENDO  
Mode



Adjustable  
Colour Temp.

MOBILE LED LIGHT



MOBILE LED LIGHT  
with battery backup



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Brightness Control	Capacitive Touch Panel
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

CLR-TA 5+5



NUMBER OF LED  
95 + 95

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,40,000 + 1,40,000 Lux ±10%



Adjustable  
Colour Temp.

**SALIENT FEATURES**  
High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Digital
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional


CLR-TA 5+4



NUMBER OF LED  
95 + 76

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
1,40,000 + 1,10,000 Lux ±10%

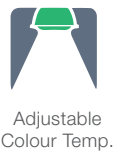

Adjustable  
Colour Temp.

**SALIENT FEATURES**  
High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Digital
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

CLR-TA 5+3



**NUMBER OF LED**  
57 + 95

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
80,000 + 1,40,000 Lux ±10%



Adjustable  
Colour Temp.

**SALIENT FEATURES**  
High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Digital
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional


CLR-TA 4+4



**NUMBER OF LED**  
76 + 76

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
1,10,000 + 1,10,000 Lux ±10%

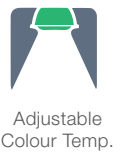

Adjustable  
Colour Temp.

**SALIENT FEATURES**  
High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Digital
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

CLR-TA 4+3



NUMBER OF LED  
57 + 76

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
80,000 + 1,10,000 Lux ±10%



Adjustable  
Colour Temp.

**SALIENT FEATURES**  
High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Digital
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional


CLR-TA 3+3



NUMBER OF LED  
57 + 57

COLOUR TEMPERATURE  
4000-5000 K

LIGHT INTENSITY  
80,000 + 80,000 Lux ±10%

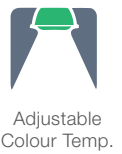

Adjustable  
Colour Temp.


**SALIENT FEATURES**  
High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm x 2
Brightness Control	Digital
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

CLR-TA 5

**SALIENT FEATURES**  
High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

- NUMBER OF LED  
95
- COLOUR TEMPERATURE  
4000-5000 K
- LIGHT INTENSITY  
1,40,000 Lux ±10%

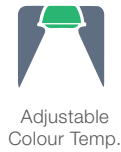

MOBILE LED LIGHT      MOBILE LED LIGHT  
with battery backup



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Brightness Control	Digital
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

CLR-TA 4

**SALIENT FEATURES**  
High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

- NUMBER OF LED  
76
- COLOUR TEMPERATURE  
4000-5000 K
- LIGHT INTENSITY  
1,10,000 Lux ±10%



MOBILE LED LIGHT      MOBILE LED LIGHT  
with battery backup



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Brightness Control	Digital
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional




CLR-TA 3

**SALIENT FEATURES**  
High quality sterilizable handle, With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

**NUMBER OF LED**  
57

**COLOUR TEMPERATURE**  
4000-5000 K

**LIGHT INTENSITY**  
80,000 Lux ±10%

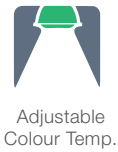

**MOBILE LED LIGHT**      **MOBILE LED LIGHT**  
with battery backup



Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Brightness Control	Digital
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

AL-HEXA

**NUMBER OF LED**  
19

**COLOUR TEMPERATURE**  
4800 K

**LIGHT INTENSITY**  
45,000 Lux ±10%



**SALIENT FEATURES**  
With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional

AL-C

**NUMBER OF LED**  
12

**COLOUR TEMPERATURE**  
4800 K

**LIGHT INTENSITY**  
40,000 Lux ±10%




**SALIENT FEATURES**  
With protective fire resistant ESG safety glass for scratch proof, fast disinfection process, protects UV and IR rays and improve light penetration. With special European design for ventilation of heat. Body is made of high quality fire resistant polymer.

Size of Light Field	12-30 cm
Colour Reduction Index	93 RA
LED Life	>50000 Hrs.
Diameter of Light	500mm
Power Supply	220V/50Hz AC
Focusing	Adjustable
Battery Backup	Optional





ISO 9001:2008 & 13485:2012

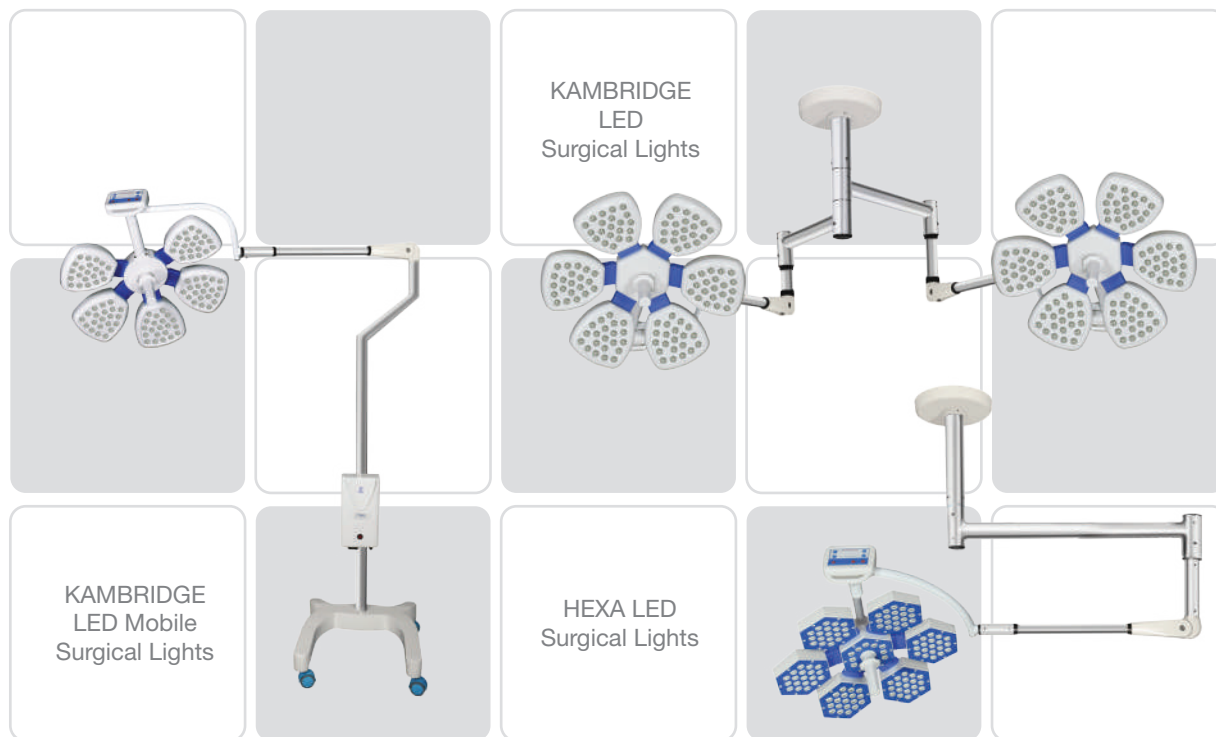

# Kay & Company

Head Office:  
25, Netaji Subhash Marg,  
New Delhi - 110002. INDIA  
Ph: +91-11-43195600 - 20 Lines  
Fax: +91-11-23275086

Works:  
B-316, Okhla Industrial Area,  
Phase - 1, New Delhi - 110020. INDIA  
Ph. / Fax : +91-11-41613432 / 33  
Email: [info@kay-co.com](mailto:info@kay-co.com), [sales@kaycoindia.com](mailto:sales@kaycoindia.com)

Representative:

For more info, please visit: [www.kaycoindia.com](http://www.kaycoindia.com)